



DEDICATED TO DISCOVERY. INSPIRED BY INNOVATION.



**HAWKINS**  
LABORATORY FUME HOODS





## OPTIONS/ACCESSORIES

<b>WORK SURFACES</b>	<b>04</b>
Epoxy Stainless Steel Type 304 Stainless Steel Type 316	
<b>BASE CABINETS - TABLE FRAMES</b>	<b>05-06</b>
General Storage Acid Storage Flammable Storage Vacuum Pump	
<b>MECHANICAL SERVICE FIXTURES</b>	<b>07-12</b>
Fixtures Pre-piping Acid Storage Vent Cabinet Vent Cup Sinks and Sinks	
<b>DILUTION TANKS</b>	<b>14</b>
2 Gallon 3 Gallon 5 Gallon	
<b>LOW AIR FLOW ALARMS</b>	<b>15-17</b>
AFA500 AFA1000 AFA4000	
<b>DISTILLATION RACKS</b>	<b>18</b>
<b>CEILING CLOSURE PANELS</b>	<b>19</b>
FINISHED BACKS SASH STOP	

## | WORK SURFACES



### Epoxy

915mm (36") wide fume hood	36-BMEP
1220mm (48") wide fume hood	48-BMEP
1525mm (60") wide fume hood	60-BMEP
1830mm (72") wide fume hood	72-BMEP
2440mm (96") wide fume hood	96-BMEP

### Stainless Steel

- Work surfaces are fabricated from Type 316 OR Type 304 stainless steel, 16 gauge, number 4 finish.
- A raised 10mm (3/8") anti-spill edge is on all four sides.
- Sinks should be fabricated from the same material and can be integral welded into the work surface.
- It shall be reinforced on the underside by 16 gauge galvanized-steel channels spaced to prevent twisting, oil-canning or buckling.
- Overall thickness shall be 32mm (1 1/4").
- Material offers excellent heat and solvent resistance and good chemical resistance to most acids.
- Stainless steel is not recommended for use with chemicals such as 18 Hydrochloric Acid, Hydrofluoric Acid, and Sulphuric Acid to 80% solution. Work surface shall be factory installed.
- Work surface is factory installed.

### Epoxy

- Work surfaces are moulded in one piece from a modified epoxy resin with a raised 10mm (3/8") anti-spill edge on all four sides and drip groove on the underside of the front edge.
- Front edge is radiused for more efficient airflow
- Overall thickness shall be 32mm (1 1/4").
- Material will offer good stain resistance and excellent chemical heat and moisture resistance.
- Sinks should be fabricated from the same material.
- Colour of top shall be black unless otherwise specified.
- Work surface is factory installed.



### Stainless Steel – type 316

915mm (36") wide fume hood	36-BMSS-6
1220mm (48") wide fume hood	48-BMSS-6
1525mm (60") wide fume hood	60-BMSS-6
1830mm (72") wide fume hood	72-BMSS-6
2440mm (96") wide fume hood	96-BMSS-6

### Stainless Steel – type 304

915mm (36") wide fume hood	36-BMSS-4
1220mm (48") wide fume hood	48-BMSS-4
1525mm (60") wide fume hood	60-BMSS-4
1183mm (72") wide fume hood	72-BMSS-4
2440mm (96") wide fume hood	96-BMSS-4

**General Storage Cabinets, Steel Construction (-S)**

- Are non-lined and have a full depth adjustable shelf.
- Cabinets are all-welded steel construction, are rigid and self supporting.
- The minimum metal gauges used in the casework fabrication shall be as follows and as recognized by North American standards:
  - o 11 gauge (3.2mm) for top and front rails
  - o 14 gauge (1.7 mm) for leveling devices
  - o 18 gauge (1.2mm) for all remaining cabinet component including cabinet frame and shelves.
- Leveling devices are installed on each corner and are accessible from the inside of the cabinet
- Cabinets are 535mm (21") deep.
- Also available 460mm deep (18").
- Also available in wood veneer construction (-W).



	880mm (34 3/4") high	825mm (32 1/2") high
	Stand-up Height	ADA Height
380mm wide (15")	No. A1-15 -S	No. ADA1-15 -S
460mm wide (18")	No. A1-18 -S	No. ADA1-18 -S
610mm wide (24")	No. A1-24 -S	No. ADA1-24 -S
760mm wide (30")	No. A2-30 -S	No. ADA2-30 -S
915mm wide (36")	No. A2-30 -S	No. ADA2-36 -S
1220mm wide (48")	No. A2-48 -S	No. ADA2-48 -S

	880mm (34 3/4") high	825mm (32 1/2") high
	Stand-up Height	ADA Height
460mm wide (18")	No. A25-18 -S	No. ADA25-18 -S
610mm wide (24")	No. A25-24 -S	No. ADA25-24 -S
760mm wide (30")	No. A25-30 -S	No. ADA25-30 -S
915mm wide (36")	No. A25-30 -S	No. ADA25-36 -S
1220mm wide (48")	No. A25-48 -S	No. ADA25-48 -S



**Acid/Corrosive Storage Cabinets, Steel Construction (-S)**

- Are designed and constructed for the storage of acids and corrosive chemicals, not flammable or combustible chemicals. Door fronts to be labeled "ACIDS" or "CORROSIVE".
- Cabinet exteriors are all-welded steel construction, are rigid and self-supporting.
- Interior lining is all welded one piece construction fabricated of 6mm (1/4") thick from white low-density polyethylene. Interior of doors are lined with the same material and have upper and lower vent slots.
- Shelf is reinforced, full depth, two position, high density polyethylene.
- The floor of the cabinet has a 25mm (1") anti-spill front edge.
- Cabinet will be provided with 40mm (1 1/2") PVC fittings and tubing for venting the cabinet into the fume hood.
- Leveling devices are installed on each corner and are accessible from the inside of the cabinet.
- Cabinets are 535mm (21") deep.
- Also available 460mm deep (18").
- Also available in wood veneer construction (-W).

## | BASIC CABINETS

	880mm (34 3/4") high	825mm (32 1/2") high
	Stand-up Height	ADA Height
460mm wide (18")	No. A23-18 -S	No. ADA23-18 -S
610mm wide (24")	No. A23-24 -S	No. ADA23-24 -S
760mm wide (30")	No. A23-30 -S	No. ADA23-30 -S
915mm wide (36")	No. A23-30 -S	No. ADA23-36 -S
1220mm wide (48")	No. A23-48 -S	No. ADA23-48 -S




### Solvent/Flammable Storage Cabinets, Steel Construction (-S)

- Designed for the storage flammable and combustible liquids and not for the storage of acids and corrosives.
- Door fronts to be labeled "FLAMMABLE- KEEP FIRE AWAY"
- Cabinets meet O.S.H.A. Standard 1910-106 (d)(3) and comply with NFPA 30 Flammable and Combustible Liquids. They are UL 1275/UL1275C UL labeled.
- Cabinets is fabricated from 1.2mm (18 gauge) steel with a baked electrostatic powder coating. The top, sides, floor and doors are double wall construction with fire-proof insulation between providing a 40mm (1 1/2") insulating air space all around.
- The floor is recessed 50mm (2") to contain spills.
- Doors are manual-closing with lever handle a three point locking devise. Self closing doors are optional
- Upper and lower air vents with spark arrestors and removable threaded covers are installed on the back of the cabinets.
- Cabinet is painted Caution Yellow.
- Cabinets are 535mm (21") deep. • Also available 460mm deep (18").
- Also available in wood veneer construction. (-W)

**NOTE:** It is not recommended to vent the cabinets into fume hood or fume hood exhaust ducting.

### Pump Cabinets, Steel Construction (-S)

- Construction is to UL962A/CSA22.2 #203 certified.
- Leveling devices are installed on each corner and are accessible from the inside of the cabinet.
- Inside of cabinet (back, sides and top) are lined with CB-300 sound deadening insulation. To protect the sound deadening insulation, a perforated painted steel lining is welded to the cabinet doors, back, sides and top.
- A 100lb capacity full extension pull out pan on the bottom of the cabinet is standard.
- One duplex electrical outlet 120V/15-20 A is located on the inside back of the cabinet, one 15A pump switch with a 4.90A thermal overload is located on the exterior top panel of the cabinet.
- A 2" (51mm) diameter vent port at the top left corner of the cabinet back and a 1 1/2" (38mm) pipe on the back right side of the cabinet top for end-user supplied piping is standard.
- Cabinets are 535mm (21") deep.
- Optional "CPU" for additional venting available
- Also available 460mm deep (18").
- Also available in wood veneer construction (-W).

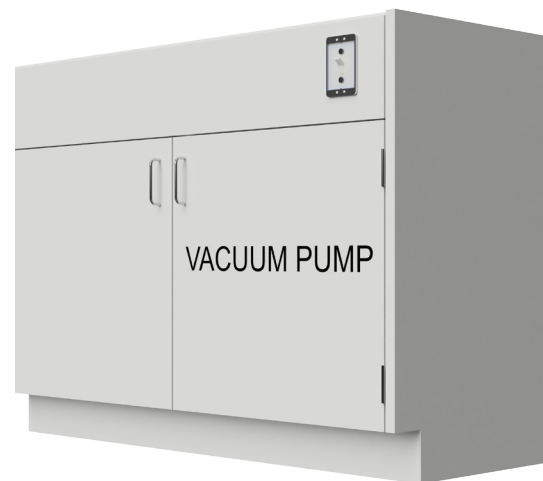

	880mm (34 3/4") high	825mm (32 1/2") high
	Stand-up Height	ADA Height
380mm wide (15")	No. A41-15 -S	No. ADA41-15 -S
460mm wide (18")	No. A41-18 -S	No. ADA41-18 -S
610mm wide (24")	No. A41-24 -S	No. ADA41-24 -S
760mm wide (30")	No. A41-30 -S	No. ADA41-30 -S
915mm wide (36")	No. A41-30 -S	No. ADA41-36 -S
1220mm wide (48")	No. A41-48 -S	No. ADA41-48 -S

**Factory installed service fixtures and fittings:**

- Mechanical service fixtures and fittings are front loading type, factory installed and piped from valve to outlet in the fume hood.
- Rod type valve assemblies are also available
- Goosenecks and nozzles are colour coded to their handles in accordance with SEFA7 - 2018.
- Fittings are CSA and CGA approved

**Valves**


- Valve Bodies: One piece construction, cast or forged brass with a minimum copper content of 57%.
- Cold Water: Ceramic compression ball valve with double sealed Uniflex connection, 1X180 degree turn open to close.
- Technical Gases (except vacuum): Brass needle valve, 3X360 degree turn open to close. Double sealed Uniflex connection.
- Vacuum: Brass needle valves, 1/4 turn open to close function. Double sealed Uniflex connection.
- Burning gases: Brass needle valve, 1/4 turn open to close with Push-In POP-UP safety handle for visual safety.

**Fixtures exposed to hood interior:**

- One piece construction, cast or forged brass with a minimum copper content of 57%.
- Finish is a chemical resisting polyester powder coating. Colour is white/ grey.

**Handles:**

- Remote control handles: Polypropylene construction four-arm handle colour coded conforming to SEFA7-2018.
- Accessible Use (ADA): One hand operated remote control handles with brass ball valve, 1/4 turn open to close. Maximum pressure to active shall not exceed 5 pounds (22.2N).

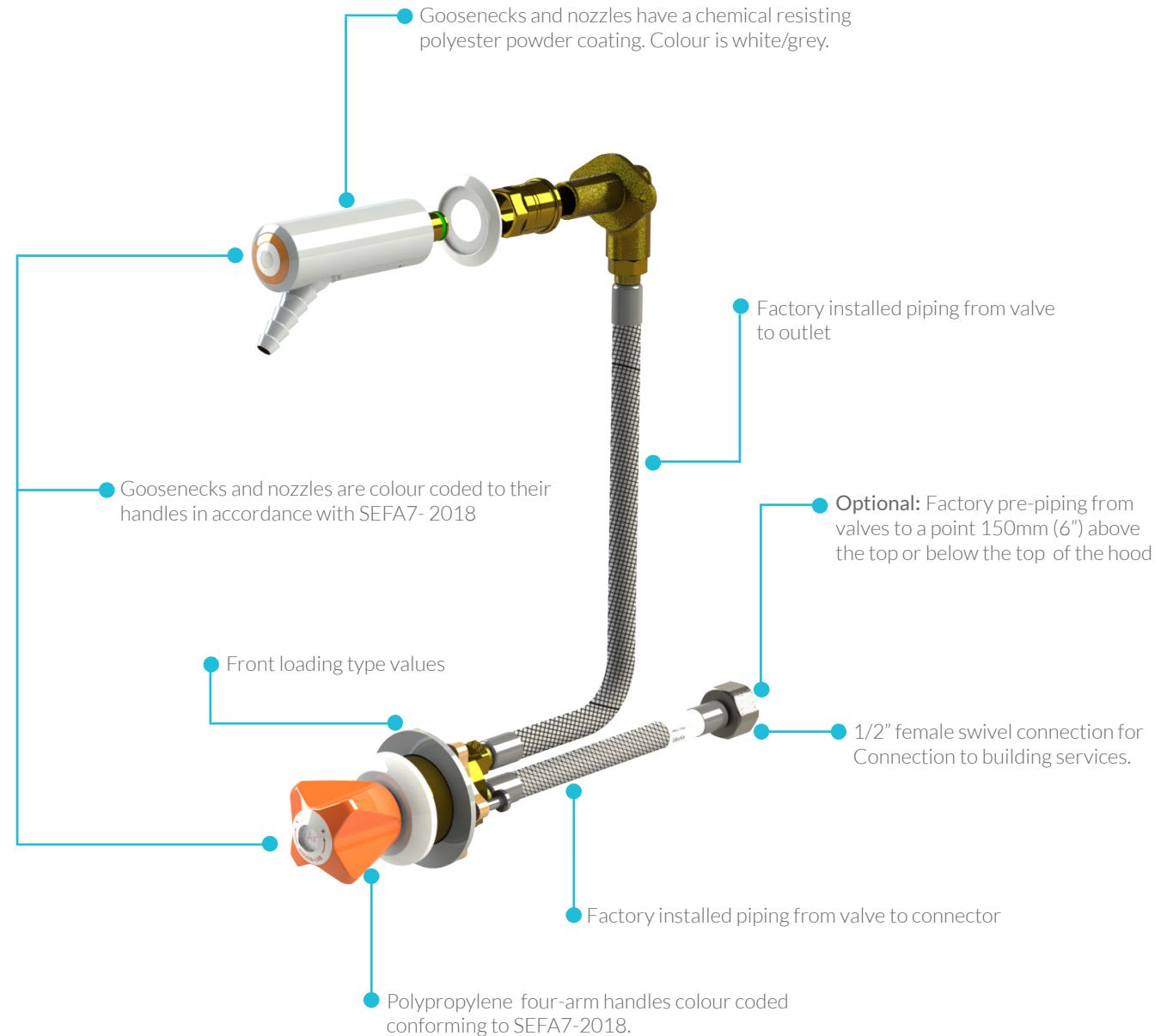

## | MECHANICAL SERVICE FIXTURES WATER



600-FL Cold water gooseneck  
602-FL Hot water gooseneck



601-FL Cold water gooseneck  
w/ vacuum breaker  
603-FL Hot water gooseneck  
w/ vacuum breaker



606-FL Cold water nozzle  
607-FL Hot water nozzle  
608-FL Cold/hot water mixing nozzle  
609-FL Steam nozzle



604-FL Cold/hot water  
mixing gooseneck



605-FL Cold/hot water mixing  
gooseneck w/ vacuum breaker



610-FL Atmospheric  
vacuum breaker





620-FL Burning gas



621-FL Air



622-FL Vacuum



623-FL Nitrogen

- 624-FL Carbon Dioxide
- 625-FL Hydrogen
- 626-FL Argon
- 627-FL Oxygen
- 628-FL Propane
- 629-FL Helium
- 630-RT Low Vacuum
- 635-FL Special gas (specify)
- 636-FL Special gas (specify)

## | MECHANICAL SERVICE FIXTURE DISTILLED/REVERSE OSMOSIS WATER

### Factory installed service fixtures and fittings:

- Mechanical service fixtures and fittings are front loading type, factory installed and piped from valve to outlet in the fume hood.
- Piping for low temperature water, 0°C to 65°C (Polyethylene and Polypropylene valves)
- Piping for high temperature water, 0°C to 90°C (Stainless steel valves and hoses)
- Goosenecks and nozzles are colour coded to their handles in accordance with SEFA7 - 2018.
- Fittings are CSA approved

### Valves

- Valve Bodies: One piece construction, cast or forged brass with a minimum copper content of 57%.
- Ceramic compression ball valve?? with double sealed Uniflex connection, 1,5X360 degree turn open to close.

### Material in contact with water:

- Polypropylene for low temperature water, 0°C to 65°C
- Type 316 stainless steel for high temperature, 0°C to 90°C

### Fixtures exposed to hood interior:


- One piece construction, cast or forged brass with a minimum copper content of 57%.
- Nozzles:
  - o Polypropylene for low temperature water, 0°C to 65°C
  - o Type 316 stainless steel for high temperature, 0°C to 90°C

### Finish:


- Chemical resisting polyester powder coating. Colour is white/ grey.

### Handles:

- Remote control handles: Polypropylene construction four-arm handle colour coded conforming to SEFA7-2018.
- Accessible Use (ADA): One hand operated remote control handles. 1/4 turn open to close. Maximum pressure to active shall not exceed 5 pounds (22.2N)..

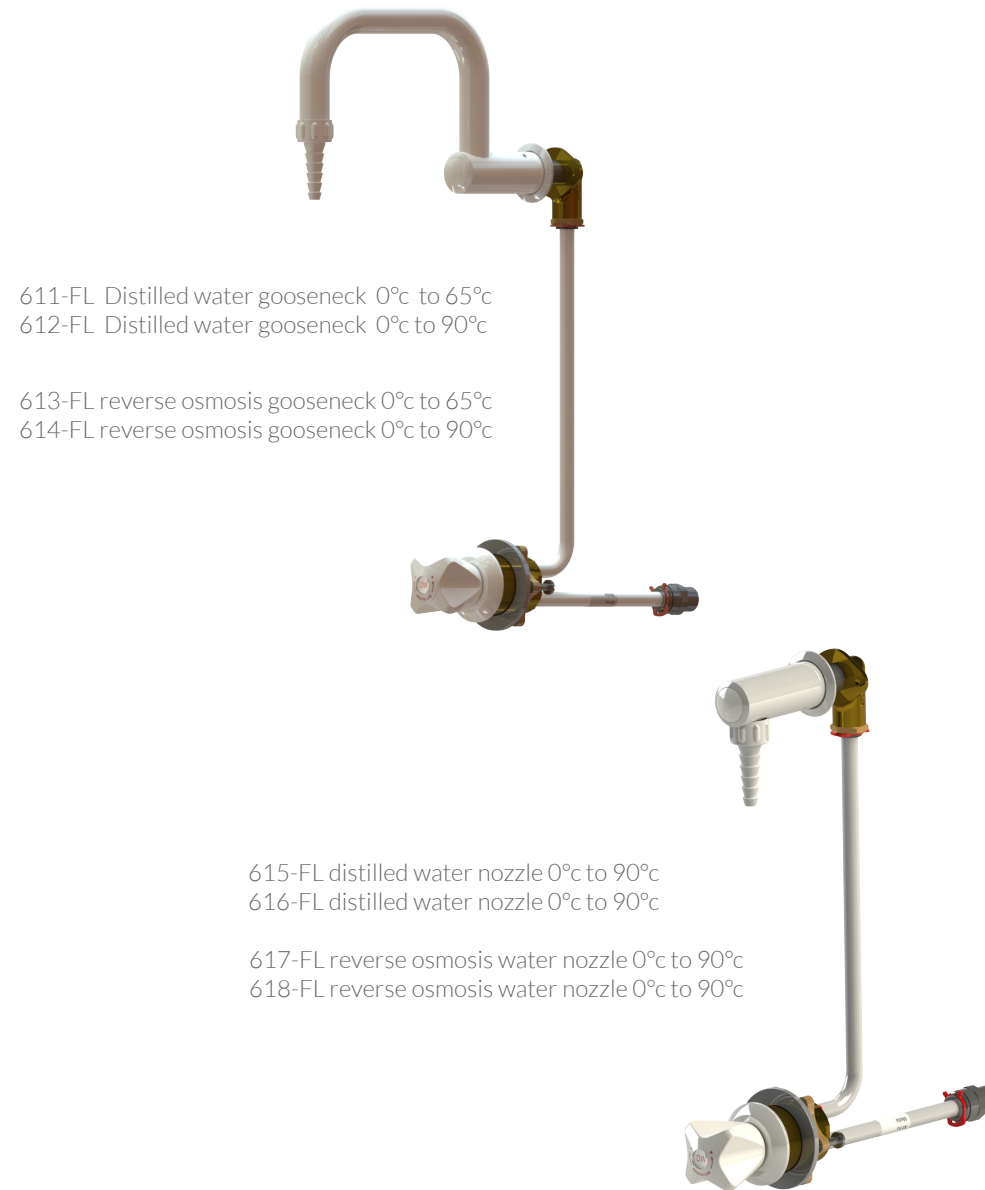

## | FACTORY PRE-PIPING | ACID STORAGE VENT KIT

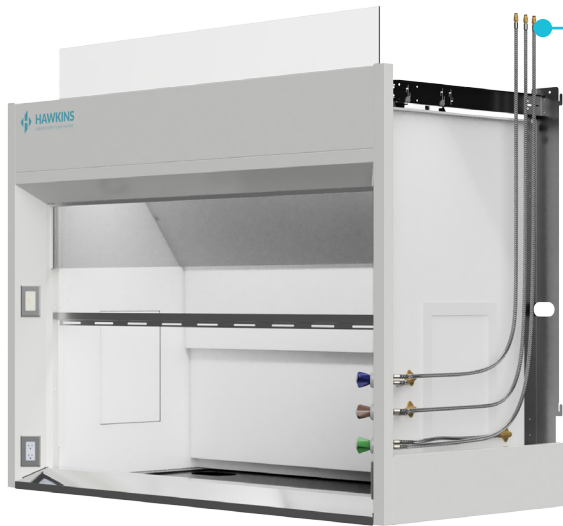

Pre-piping 150 (6")  
Above Hood

### Factory Pre-Piping - Optional

- Pre-piping from valves to a point 150mm (6") above OR below the hood superstructure complete with 1/2" female nut connector for service connection.

### Material:

- Water: SPX (Soft PEX-A) hose with stainless steel braiding. Outside diameter is 0.47" with a minimum inner bending radius of 1 37/64".
- Technical gases: SPX (Soft PEX-A) with (describe) stainless steel braiding. Outside diameter is 0.47" with a minimum inner bending radius of 1 37/64".
- Burning Gases: Corrugated type 316 Stainless steel hose. Outside diameter is 0.47" with a minimum inner bending radius of 0.33".

### Pressure Ratings:

- Fittings are to be constructed to operate with the following maximum working pressure without leak or failure.
  - Water Fittings: 145 PSI
  - Non-Burning Gas: 145 PSI - 231 PSI
  - Burning Gases: 100 PSI
  - Special Water Fittings: 145 PSI
  - Oxygen Fittings: 145 PSI - 232 PSI



Pre-piping 150 (6")  
Below Hood



### Acid Storage Cabinet Vent Kit

- Designed to vent the acid storage base cabinet into the rear upper baffle area of the fumehood without perforating the work surface.
- Assembly consists of 40mm (1 1/2") diameter PVC tubing and PVC fittings.
- Tubing is connected to the back of the cabinet and the rear baffle area of the hood with PVC elbows and flanges.
- PVC tubing runs between the inside and outside wall of the fume hood superstructure.

## | EPOXY SINKS AND POLYETHYLENE CUP SINKS

### EPOXY SINKS AND CUP SINKS

- Moulded seamless epoxy construction. Standard colour is black. Grey and slate grey are optional colours.
- Sinks are Drop-in type with moulded rim. Can be flush mounted in the counter top or with raised rim 1/4" above the countertop surface.
- Waste fitting is polyethylene with 1 1/2" diameter threaded tail piece and strainer. Overflow tube is optional.




White, grey and slate grey are available to order with no up-charge.

Model	Internal Length	Internal Width	Internal Depth	External Height	Flange Width	Internal Thickness
	A	B	C	D	E	F
#3	12	8	8	9	0.8	0.3
#5A**	14	10	6	7	0.8	0.3
#10**	16	8	7	8	0.8	0.3
#15**	16	12	8	9	0.8	0.3
#25A	18	15	5	6	1.1	0.3
#25	18	15	8	9	1.1	0.3
#30**	18	15	11	12	0.8	0.3
#33**	21	17	10	11	0.8	0.3
#50	24	16	8	9	0.8	0.3
#55A	25	15	5	6	0.8	0.3
#55**	25	15	10	11	0.8	0.3
#59	28	15	12	13	0.8	0.3
#70	24	16	16	17	0.8	0.3

Available in black.

Model	Internal Length	Internal Width	Internal Depth	External Height	Flange Width	
S497	9	3	5.8	6	0.63	Deck Mount
S499	6	3	4.5	5	0.63	Deck Mount
S500	5.4" Diameter	5.4" Diameter	3.7	4	0.63	Deck Mount
S502	6	3	3.7	4	0.63	Wall Mount

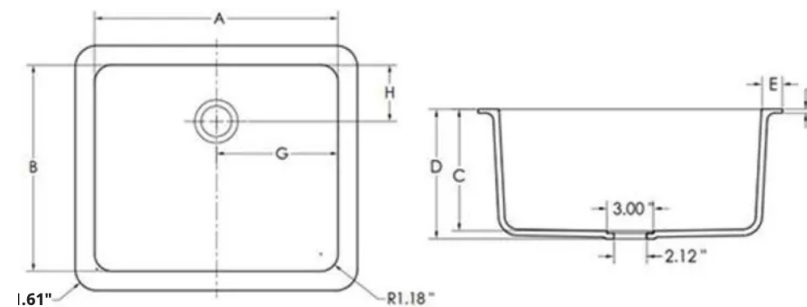

- Stamped seamless construction, type 316, 20 gauge. Number 4 satin finish on exposed surfaces.
- All corners are radiused for ease of cleaning.
- Self rimming with integral rim. Can be integrally welded into stainless steel work surface.
- Undercoated for sound deadening and reduced condensation
- Waste fitting is a 3 1/2" crumb cup stainer with stopper and 1 1/2" diameter stainless steel tailpiece




**INSIDE DIMENSION**

NO.SS1	14½" x 9½" x 7"
NO.SS3	16" x 14" x 10"
NO.SS4	18" x 16" x 7"
NO.SS5	18" x 16" x 10"



**INSIDE DIMENSION**

NO.SS30	6" x 3" Oval
NO.SS30	6" Diameter

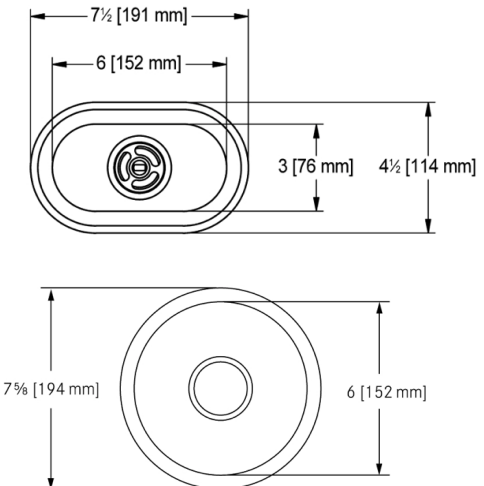


- Stamped seamless construction, type 316, 20 gauge. Number 4 satin finish on exposed surfaces.
- All corners are radiused for ease of cleaning.
- Self rimming with integral rim. Can be integrally welded into stainless steel work surface.
- Undercoated for sound deadening and reduced condensation  
Stainer is built in with 1 1/2" diameter stainless steel tailpiece

## | DILUTION TANKS

### DILUTION TANK/TRAP - 2 GALLONS

- Design allows for easy fitting in tight spaces.
- Standard connections are 1-1/2" female threaded inlet and 1-1/2" male threaded outlet.
- A 3" cleanout is mounted on the tank top to facilitate inspection and maintenance.
- Single or dual inlets can be provided. The standard tank configuration is shown.
- May be used as sedimentation traps.

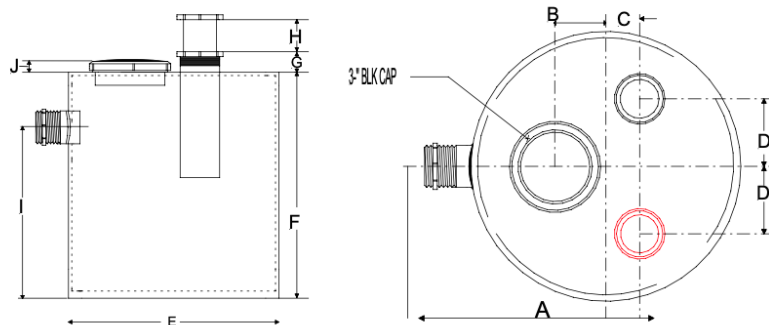


WHEN TANK IS ORDERED WITH 1 INLET  
THE STD LOCATION IS MARKED IN RED.

COLOR	PART #	A	B	C	D	E	F	G	H	I	J	INLET	# OF INLET	OUT-LET
BLACK	W59155-200	7-5/16"	2-1/2"	3/16"	2-1/16"	9"	10-5/32"	1-15/16"	0-10"	6-13/16"	1-5/8"	1-1/2"	2	2"
BLACK	W59155-201	7-5/16"	2-1/2"	3/16"	2-1/16"	9"	10-5/32"	1-15/16"	0-10"	6-13/16"	1-5/8"	1-1/2"	1	2"
BLACK	W5915-202	7-5/16"	2-1/2"	3/16"	2-1/16"	9"	10-5/32"	1-15/16"	0-10"	6-13/16"	1-5/8"	2"	1	2"
BLACK	W59155-203	7-5/16"	2-1/2"	3/16"	2-1/16"	9"	10-5/32"	1-15/16"	0-10"	6-13/16"	1-5/8"	2"	2	2"

ALSO AVAILABLE IN TRANSLUCENT WHITE W59156


### BOOKSHELF/BOILER TANK

Confined space tanks all dimensions are in inches

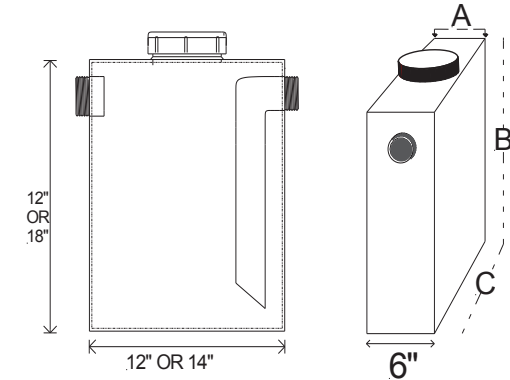

ONE PIECE ROTATIONALLY MOLDED  
1 1/2" OR 2" MALE/FEMALE  
OR PLAIN END INLET AND OUTLET  
3" TOP CLEANOUT CAP

A	B	C	D= GALLONS	PART NUMBER
6"	12"	12"	3	W121206-001
6"	18"	14"	5	W181406-001

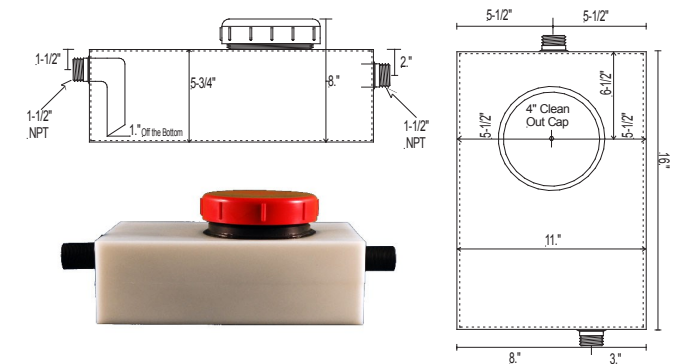



## AFA 500

Model AFA500 factory installed low airflow monitor and alarm system.

- Built in sidewall airflow sensor for measuring and monitoring face velocity.
- Audible alarm with mute button.
- Visual alarm with Green LED = Safe and Red LED light = Alarm.
- Audible alarm can be silenced but Red LED will stay active until alarm condition is corrected.
- Single Alarm Point Capture pushbutton calibration.
- Relay input for Night Setback to mute audible alarm.
- Alarm is low voltage and is supplied with a transformer.

## | LOW AIRFLOW ALARM MONITOR



### AFA 1000

Model AFA1000 factory installed low airflow monitor and alarm system

- Built in sidewall airflow sensor for measuring and monitoring face velocity
- Back-lit LCD display.
- Digital face velocity display in fpm or m/s.
- Visual and audible alarms.
- Indicator display: Green LED = Safe, Amber LED = Caution, Red LED = Alarm
- Alarm Indiction: Red graphic with audible alarm.
- Audible alarm can be silenced but Red graphic will stay active until alarm condition is corrected.
- Two point alarm capture push button calibration.
- BACnet and Modbus on board available with optional comms adaptor.
- Alarm is low voltage and is supplied with a transformer.





## AFA 4000

### Model AFA4000/1 factory installed low airflow monitor and alarm system

- Built in sidewall airflow sensor for measuring and monitoring face velocity
- Full-colour 3.5" LCD display.
- Digital face velocity display in fpm or m/s.
- Visual and audible alarms.
- Graphic display: Green LED = Safe, Amber LED = Caution, Red LED = Alarm
- Alarm Indiction: Red graphic with audible alarm.
- Audible alarm can be silenced but Red graphic will stay active until alarm condition is corrected.
- Two point alarm capture pushbutton calibration.
- BACnet and Modbus on board available with optional comms adaptor.
- Alarm is low voltage and is supplied with a transformer.
- Compatible with TEL Variable Air Volume Auto Sash Return System

## | DISTILLATION RACKING

### Assemblies

- Nominal 300mm X 300mm (12" X 12") grid pattern
- Adjustable by users as apparatus requirement change
- Available in three materials to better suit user needs
- Available for bench mount and floor mount hoods
- Custom sizes to customer requirements

### Stainless Steel, type 316, satin finish (S)

- ½" diameter stainless steel rounded edge type rods
- Open end style stainless steel clamps with set screw
- Stainless steel flanges with set screw


### Aluminum, satin finish (A)

- ½" diameter aluminum rounded edge type rods
- Open end style aluminum clamps with set screw
- Aluminum flanges with set screw


### Plastic (P)

- ½" diameter white fiberglass rounded edge type rods
- Open end style black nylon clamps with set screw
- Black nylon flanges with set screw


## BENCH MOUNT

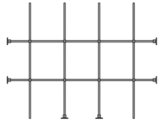


BMDR-48(A)  
BMDR-48(S)  
BMDR-48(P)

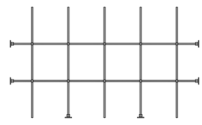

BMDR-60(A)  
BMDR-60(S)  
BMDR-60(P)

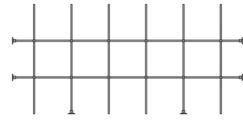


BMDR-72(A)  
BMDR-72(S)  
BMDR-72(P)

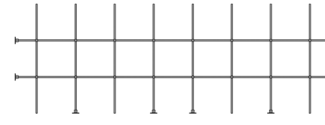


BMDR-96(A)  
BMDR-96(S)  
BMDR-96(P)


## FLOOR MOUNT

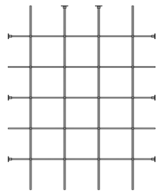


FMDR-48(A)  
FMDR-48(S)  
FMDR-48(P)

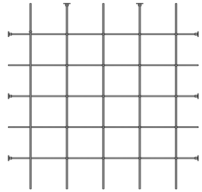

FMDR-60(A)  
FMDR-60(S)  
FMDR-60(P)

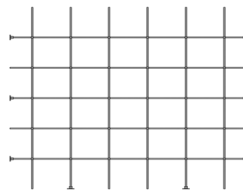

FMDR-72(A)  
FMDR-72(S)  
FMDR-72(P)

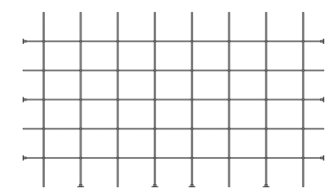


FMDR-96(A)  
FMDR-96(S)  
FMDR-96(P)





#### Ceiling Closure Panels

- Ceiling Closure Panels are intended to enclose the top of the fume hood to the underside of the ceiling, so as to hide the exhaust ducting from view.
- Fabricated from the same material and gauge as the fume hood exterior.
- Panels are a set-back design to allow for proper bypass air operation and to allow for the removal of the fume hood front and side panels.
- Maximum height is 610mm (24").
- Colour to match fume hood exterior.

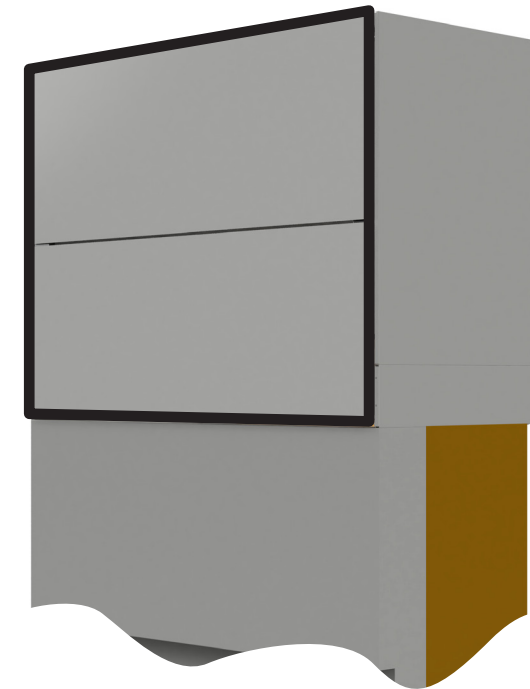


#### Finished Back Panels

- Finished Back Panels are intended to enclose the back of the fume hood when the hood is sitting on an island or peninsula cabinet assembly.
- Panels are a two-piece construction and are removable.
- Colour to match fume hood exterior.

#### Sash Stops

- Sash stops can be provided at the specified Sash Operating Height.
- They have a manual override when lowering the sash below the Sash Operating Height and an automatic reset when the sash is raised above the Sash Operating Height.

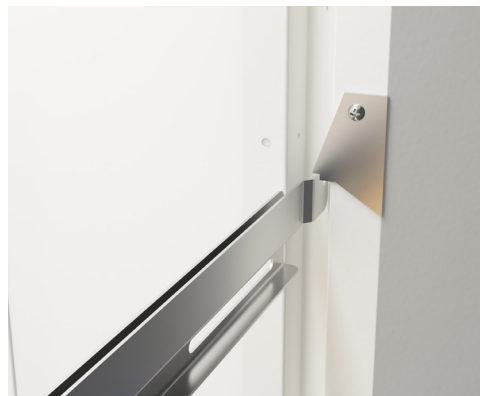

YOUR RESEARCH DRIVES OUR INNOVATION.





Unit 107 - 19298 21st Ave Surrey, BC V3Z 3M3

1-800-661-4454